

# NIDO™ Pedicle Screw System

The NIDO™ Pedicle Screw System is a comprehensive modular experience designed for a midline approach. The system features a broad spectrum of ultra-low-profile implants with a focus on surgical options for multi-level fusions.

*NIDO™ is a trademark of Kalitec Medical.*



# SYSTEM HIGHLIGHTS

*Ultra-Low-Profile Tulips  
Featuring Only  
8mm Run on Rod*


*Self-Tapping Screw Diameters  
4.75, 5.5, 6.5, 7.5, 8.5, 9.5, and 10.5mm*

*Fenestrated Screws Available Upon Request*



*Screw Lengths  
25-115mm  
(Varies by Diameter)*

*Dual-Lead Cortical  
Cancellous Thread  
Profiles*

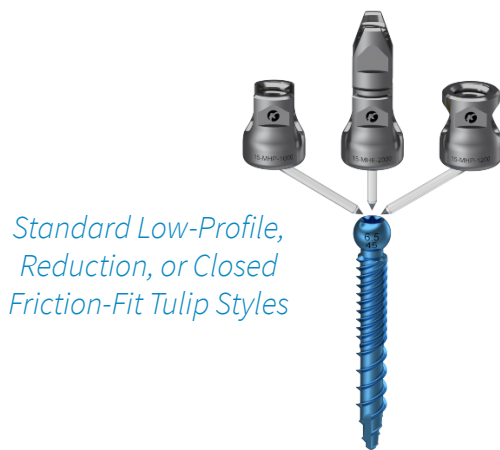

*Standard Low-Profile,  
Reduction, or Closed  
Friction-Fit Tulip Styles*

*Self-Tapping Tapered Tip for Ease of Insertion*



*Negative Thread  
Angle Prevents  
Possible Tulip Splay*

## Domino Rod-to-Rod Connectors

*Featuring Low-Profile  
Axial, Open-Open, Open-  
Closed, and Closed-Closed  
Geometries*



## Tulip Rod-to-Rod Connectors

*Modular Parallel  
Tulips Attach to  
Any Size NIDO  
Screw for Simplistic  
Implementation*



*Adjustable Parallel  
Tulips Provide Sagittal  
and Transverse  
Degrees of Freedom*

